



IBT BioServices  
20358 Seneca Meadows Parkway  
Germantown, MD 20876  
240-454-8934 x229  
Services@integratedbiotherapeutics.com

## Chikungunya 181/25 Antibody

**Catalog #:** 01-0008

**Lot #:** CH200909IBT

**Immunogen:** Sucrose purified virus from infected supernatant; CHIKV 181/25

**Description:** Rabbit polyclonal for CHIKV strain 181/25

**Raised in:** Rabbits

**Clonality:** Polyclonal

**Relevance:** The 181/25 vaccine strain represents a surrogate strain for the Category C (NIAID) virus Chikungunya.

**Recommended Dilutions:**


ELISA: Assay-dependent dilution

WB: Assay-dependent dilution; if probing against purified infected cell supernatant, cross-adsorption recommended.

IF: No Data

**Storage:** 2-3 weeks +4°C, -20°C long term

**Cross Reactivity:** Does not react with Vaccinia NYCBOH, Junin Candid #1, or Rift Valley Fever MP12. Some cross-reactivity observed with fetal bovine serum, Vero cells, and Venezuelan Equine Encephalitis strain TC-83.

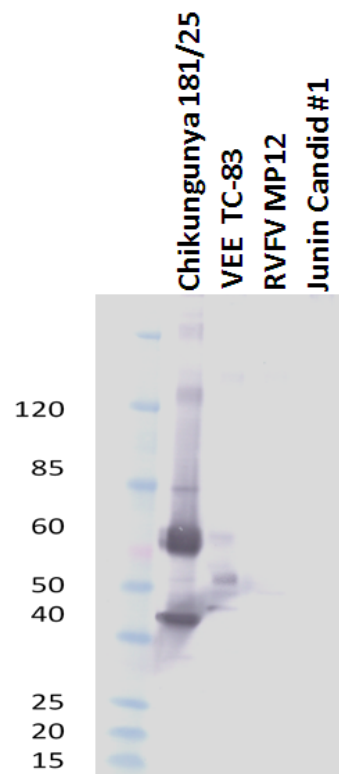

-Serum tested against sucrose purified infected cell-supernatant  
-Serum diluted 1:10,000 and incubated with 2% FBS for 2hr  
Prior to blotting

**Intended for research use only, not for human, therapeutic, or diagnostic applications.**